


# C-RCC

Root Canal Irrigation and Suction System



## Why we need to irrigate?

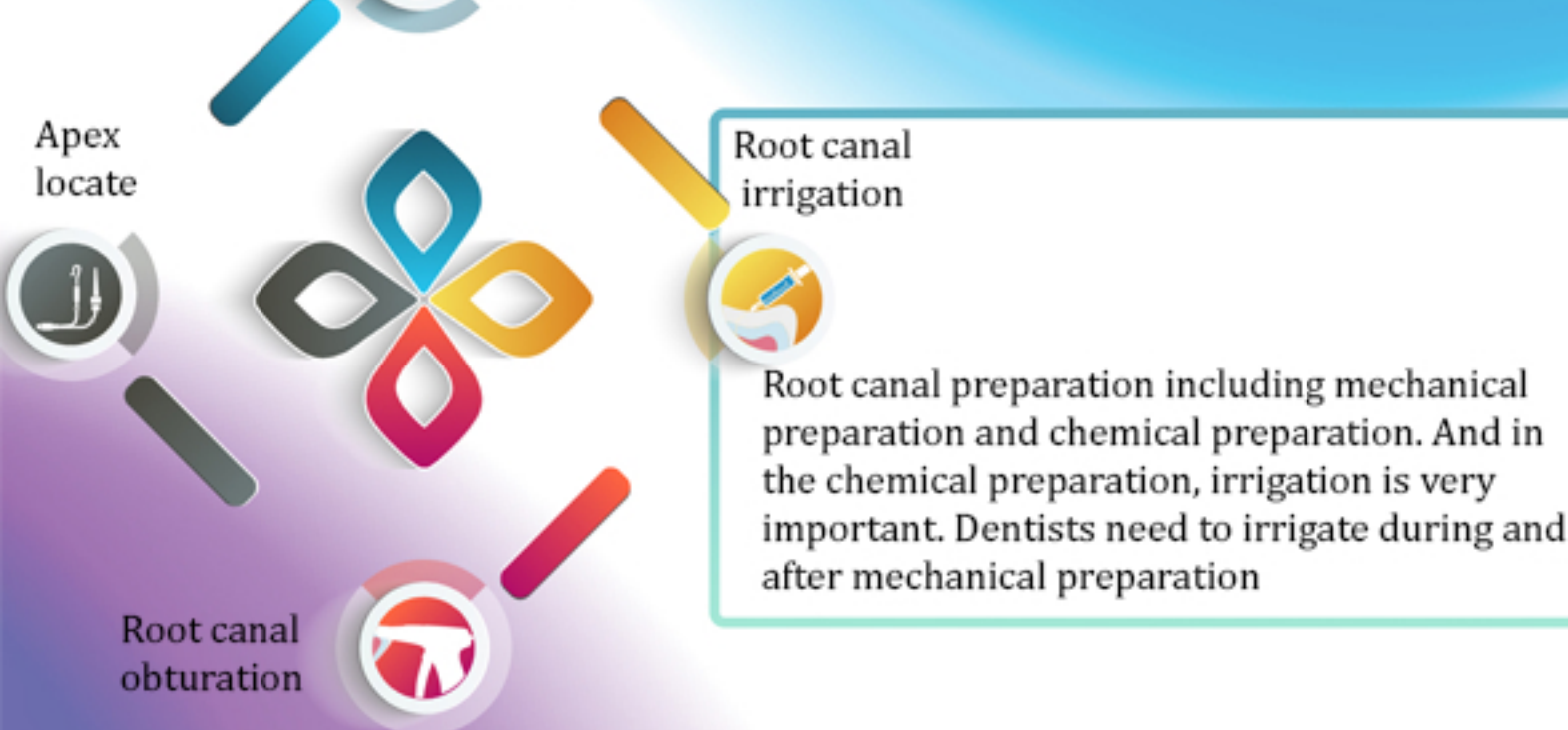

## Innovative look

Beauty & Fashion



## Multi-working function


Irrigation and suction working together.



Make your treatment more convenient


## With heating function

- Improve the disinfection effect by 400%, making the treatment safer, faster and more efficient.
- Can turn off the heating function Which can be used with normal saline and liquid EDTA

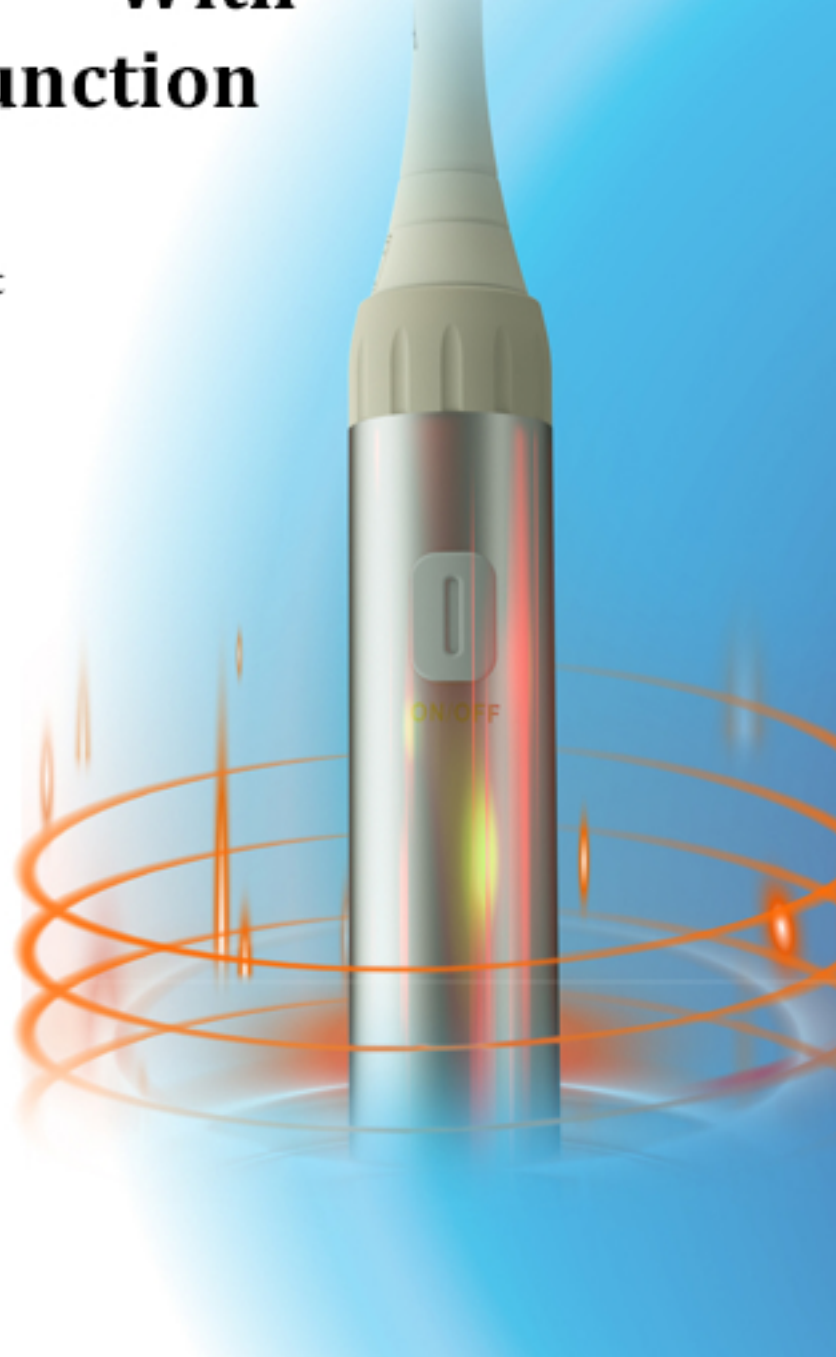

## Negative pressure


Avoid increasing apical pressure and reduce the risk of apical overflow, safer

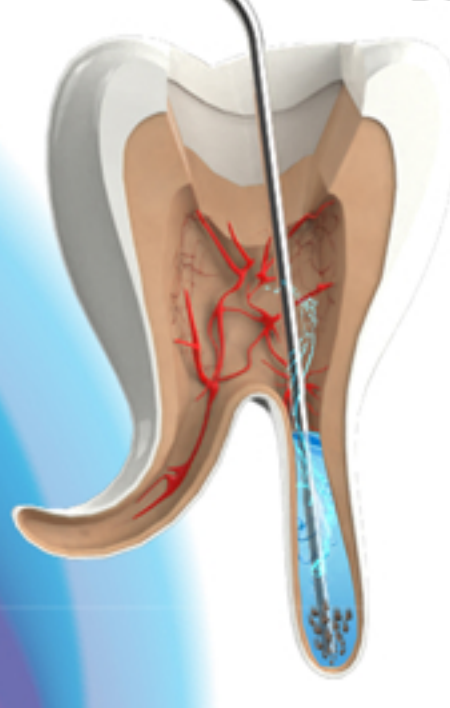

## Multi-function

Meet all your requirement

- 1.Irrigate only
- 2.Suck only
- 3.Suck while irrigating
- 4.Cleaning mode
- 5.Acclerate mode
- 6.Working time adjustable

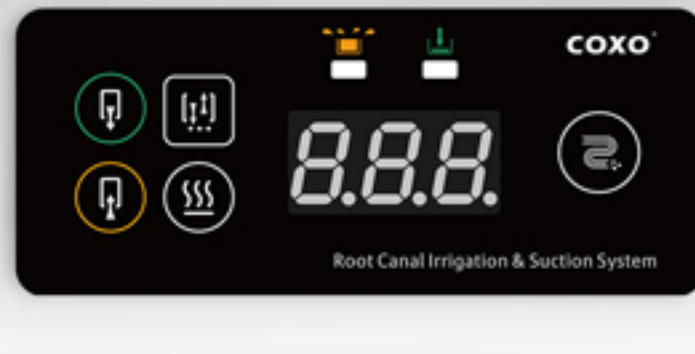

## Special tip design

### Tip A

For removing larger debris, Use after each time mechanical preparation Shorter needle is bigger, size 21G, for suction Longer needle is smaller, size 27G, for irrigation

### Tip B

For removing smaller debris Use after finishing mechanical preparations Shorter needle size 27G, for irrigation Longer needle size 27G, for suction

## Humanized design

**Different Bottle Design**  
Make sure you will not take the wrong one



**Big Volume Design**  
Enough for 15 teeth meet all your need

Adaptor working voltage:	AC110-220V 50/60Hz
Working time:	30s/60s/90s/120s
Bottle volume:	irrigation:100ml/waste :200ml
Irrigate volume:	>10ml/min
Suction volume:	>3ml/min
Volume(cm):	34x26x10
Weight(kg)/PCS:	1.65

