

Good Doctors  
Wireless Solution

Caries  
Detector

Wireless  
Intraoral Camera

Wired  
Intraoral Camera

Apex  
Locator

Curing  
Light

IC-AD

IC-WHCD100  
IC-DCDU75

IC-WHT60  
IC-DCAM70  
IC-EZAM80

IC-WHMA90 IC-EZAU81  
IC-WHW61 IC-DCCU74  
IC-DCAU71 IC-EZBU82


AL-DFA20  
AL-DF10

CL-DP40 CL-DC20  
CL-DP41  
CL-DL10



### Zoom function

Ability to enlarge the image with zoom feature to support user convenience.

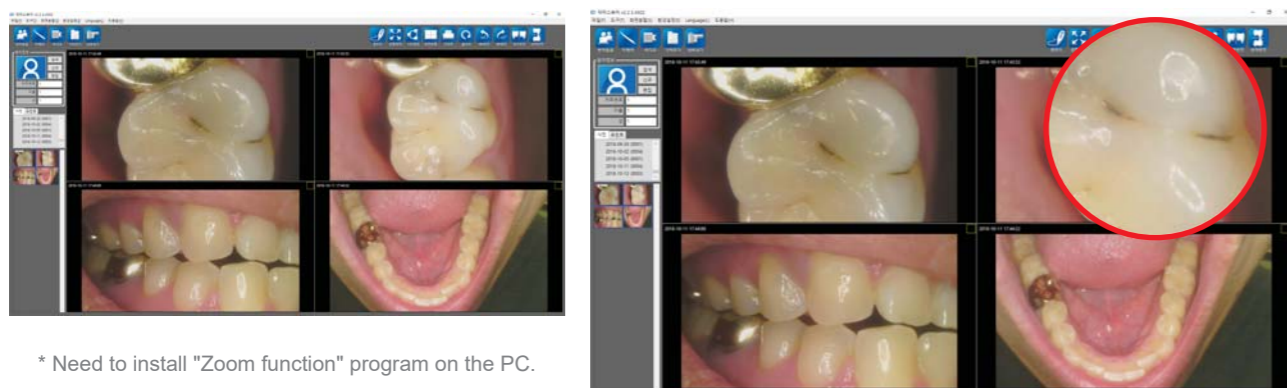

\* Need to install "Zoom function" program on the PC.



### Software compatibility

Able to be used with almost any dental imaging software in the market.



## 1 Specification

Model	IC-WHMA90	
Camera	Operation voltage	5V 600mA
	Measurement	Approx. 128 x 21.5 x 32 mm
	Weight	Approx. 61g
	Resolution	FHD (1920x1080)
	Sensor	5M Sensor CMOS / Liquid
	Focusing	Auto Focus
	Focus range	Macro to Infinity
USB cable	cable length	Approx. 3 M (9.8ft)

## 2 Recommended PC system

PC system:	3.3GHz CPU, 4GB RAM
OS:	Windows 7
USB port:	USB 2.0



Good Drs



Whicam M

Macro to Infinity

Autofocus Dental Camera with Mouse



GDBR-IC-WHMA90(Brochure for Whicam M)(Eng)(2)



www.gooddrs.com  
E-mail. info@gooddrs.com



Good Doctors Co., Ltd.  
#208, B-dong, 283 Bupyeong-daero, Bupyeong-gu (Woolim Lions Valley,  
Cheongchoen-dong), Incheon, 21315, Republic of Korea  
Tel: +82 32 424 6325 Fax: +82 32 424 6326



# Features of Whicam M

## Image quality

- Macro to infinity**  
Camera is able to focus from macro mode (close) to infinity that includes a full arch and face shot.
- Auto-focus**  
Camera is able to focus in any position not only inside of oral cavity but also full arch and face shot.
- FHD image support**  
The device shows FHD (1080p) in live image supported by a 5M sensor.
- Color impressions**  
Adjust the color impression on camera to fit different monitors or screen.
- Black&White function**  
Capture and save X-ray film image to PC



## Ergonomic design

- Dental handpiece design**  
Most familiar design for dentists.
- Metal body case**  
Metal body case increases product durability.
- Slim head size**  
Slim head design makes it easier to approach the oral cavity with camera.
- Retraction head design**  
Retract patient's tongue and cheek to get better image.
- Efficient brightness**  
12 pcs LEDs makes for efficient brightness.
- Dual capture buttons**  
Capture buttons are placed on the top and bottom side of camera to make for an easy capture in any position
- LED indicators**  
LED indicators show the device status and mode.



## Convenience

- Mouse function**  
Built-in mouse function for user convenience.
- Zoom function**  
Ability to enlarge the image with zoom feature to support user convenience.
- LED brightness adjust**  
Ability to adjust the LED brightness.
- Support twain driver**  
Support twain driver allows for use in almost any dental imaging software.
- Software compatibility**  
Able to be used with almost any dental imaging software in the market.
- Detachable USB cable**  
Detachable USB cable makes it easy to move the camera to other dental chairs.




### Macro to Infinity

Camera is able to focus from macro mode (close) to infinity that includes a full arch and face shot.



### FHD image support

The device shows FHD (1080p) in live image supported by a 5M sensor.

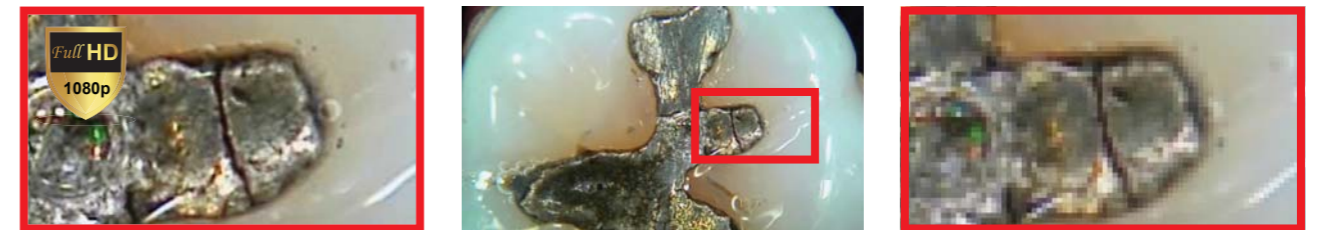

Full HD


Original image

HD



### Color impressions

Adjust the color impression on camera to fit different monitors or screen.

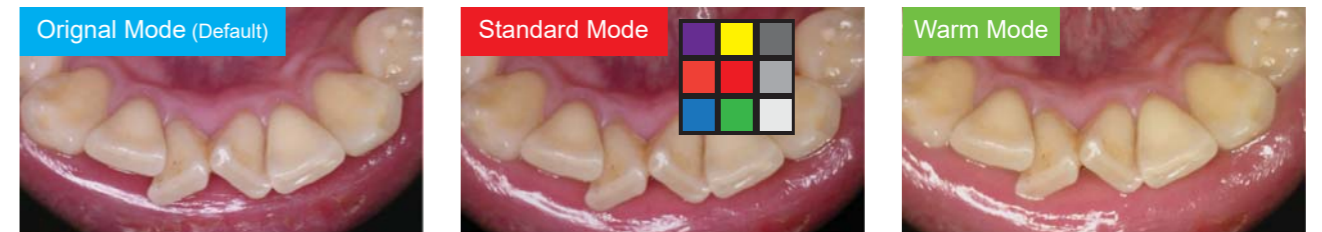

Same color chart image with xrite

