

# Portable Dust collector

Portable  
Dust  
Collector

Outstanding  
Suction  
Force

Wireless  
Sensor  
Technology



Portable Dust Collector for Dental Clinic and Lab

# Portable Dust collector




## Specification

Size	220(W) x 230(H) x 260(D)	Zabara	50Ø, max 500mm
DC Adapter	19V5A, 95W	Sensor	Apply PIR sensor
Battery capacity	4S/1P(Li-ion) 2900mA 43W	Filter	PM 2.5
Charging time	3 hour	LED (Indicator)	5,500K
FAN wind pressure	1,280Pa	Botton	Touch sensor
Voltage	AC100~240v, 50~60Hz	Output	DC 12V


### LED Lamp close to sun light at 5,500K

With a near natural light LED lamp of 5,500K, there is little eye strain even when used for long periods of time.

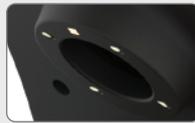

### PIR sensor prevents Cross infection

The use of the Human Body Detection Sensor (PIR) reduces button touch ans is excellent for preventing cross-infection.

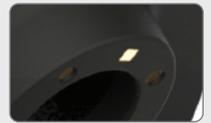

### Easy to detach with a magnet flexoible bellows (Max 500mm)

A magnet provides easy removal and a maximum length of 500mm increases the bag capacity.



### Apply ultra fine dust filter PM 2.5

We applied PM 2.5 filter to remove 99% of ultra fine dust.



Tel. +82-1899-2804 | Fax. +82-53-583-2806  
99, Seongseoseo-Ro, Dalseo-Gu, Daegu, Korea



Tel. +1-323-677-4363~5 | Fax. +1-323-677-4366  
6 Centerpointe Drive, Suite 600 La Palma CA 90623



Tel. +021-5111-3828 | Fax. +021-5111-3828  
上海市长宁区中山西路933号2205室